

DR Double Pump

Ref. No. ARTH-PUMP-A



The DR Double Pump is an integrated fluid management system (inflows and outflows).

The system automatically regulates irrigation relative to aspiration, thereby ensuring a perfusion rate and precisely controlling pressure.

Features

- Automatic regulation of inflow and outflow: Monitors perfusion and pressure
- Integrated shaver interface for most shaver brands
- Increased responsiveness for optimal convenience during use
- Pressure readings in mmHg or LH2O
- Familiar user interface
- The arthroscopic tubing sets are designed for the irrigation of joint cavities with saline solution and for aspiration to remove debris during surgical procedures.

Very user-friendly and simple to operate.

Technical parameters

Pressure range	20 – 280 mmHg (10 – 120 LH2O)
Tubing perfusion rate	At rest: 100 ml / min \pm 50 ml / min Min: 200 ml/min \pm 50 ml/min Low: 300 ml/min \pm 75 ml/min Medium: 500 ml/min \pm 75 ml/min High: 600 ml/min \pm 75 ml/min
Shaver perfusion rate	Min: 200 ml/min \pm 50 ml/min Low: 400 ml/min \pm 75 ml/min Medium: 600 ml/min \pm 75 ml/min High: 800 ml/min \pm 100 ml/min
Technical data (electrical)	Line voltage: 100 - 240 V AC Frequency: 50-60 Hz Max. power consumption 500 VA
Dimensions (W x H x D)	434 x 202 x 377 mm
Weight	11 kg

Technical parameters

DR Single-Day Tubing Set


Ref. No. AV4503

Single-day tubing: An irrigation cannula that conducts a salt solution from a bag into connector tubing.

Thanks to a one-way valve, the single-day tubing set can be re-used with patient tubing or connector tubing, which keeps your costs low.

- A single irrigation cannula for the entire workday
- Your solution bag use reduced
- Time saved in setting up new tubing

Sterile packaging, box of 30

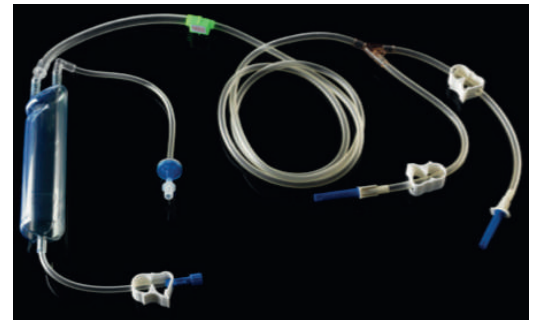

DR Connector Tubing with One-Way Valve

Ref. No. AV1102

Connector tubing with one-way valve: An irrigation cannula to be connected with the DR Single-Day Tubing Set.

If only irrigation is required, the connector cannula is used in combination with the DR Single-Day Tubing Set.

Sterile packaging, box of 50



DR Patient Tubing with One-Way Valve

Ref. No. AV4509

DR patient tubing with one-way valve: Used together with the DR Single-Day Tubing Set when irrigation and outflow control are required.

DR Patient Tubing includes:

- Connector tubing with one-way valve
- Outflow tubing to be connected with the aspiration cannula or the DR handpiece to dispose of blood and debris in a waste container.

Sterile packaging, box of 30



DR 2-Way Foot Pedal

Ref. No. ARTH-PED2-A

Intended for use in combination with shaver interface clip or integrated shaver interface.

DR 4-Way Foot Pedal

Ref. No. ARTH-PED4-A

Intended for use in combination with foot pedal interface.



Legal Manufacture: Hemodia SAS / 85 rue du Chêne vert / 31670 Labège / France / www.hemodia.com

Version001_EN