

DR Universal Single Pump

Ref. No. 200.000.034



The DR Universal Single Pump was designed to supply fluids to a procedure site at a specific pressure.


The device is for arthroscopic, hysteroscopic, and urologic operations.

Operated via control unit or DR Universal Single Pump remote control.

For use with disposable and single-day tubing.

Functions

- Manufactured to dispense liquid to a procedure site
- Settings can be saved
- Achieves efficient rinsing at a procedure site
- Acoustic and visual signals ensure that settings are correct
- Can be operated using DR Universal Single Pump remote control
- Large and clear display

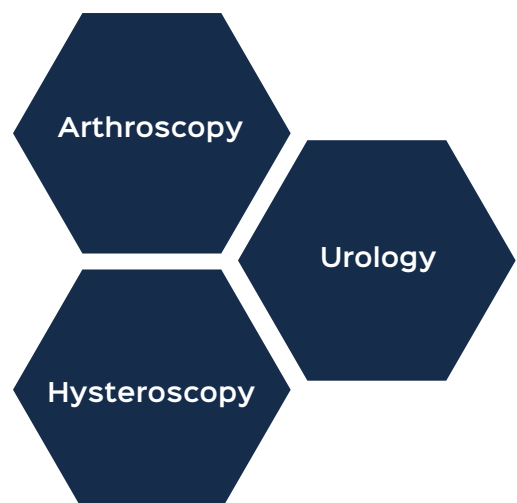

Technical parameters

DR Universal Single Pump

Ref. No. 200.000.034

Voltage	100-240 V AC
Current frequency	50 / 60Hz
Housing safety standard	IP X1
Ambient temperature	When operating: +10°C bis +40°C When in storage or during transport: -20° C bis +45°C
Dimensions (W x H x D)	330 x 120 x 306 mm
Weight	7,5 kg



DR Tubing for the DR Universal Single Pump

DR tubing sets for the DR Universal Single Pump. Designed to make it simple for medical personnel to prepare pumps for an operation. Two options are available: DR disposable tubing; DR single-day and DR patient tubing set, the more economic solution.

DR disposable tubing set

Ref. No. 200.000.044

Installed in direct contact with the pump control unit sensor.
Use of the DR Disposable Tubing Set mandates a change after every use.

Sterile packaging, box of 10



DR tubing for economical use, consisting of

DR single-day tubing

Ref. No. 200.000.036

Installed in direct contact with the pump control.

Sterile packaging, box of 10



DR patient tubing

Ref. No. 200.000.037

Installed in direct contact with DR single-day tubing and with an endoscope port.

Sterile packaging, box of 50



DR universal single pump remote control

Ref. No. 200.000.035

The DR Universal Single Pump remote control:
For the DR Universal Single Pump.

This makes it possible to control liquid flow and pressure without being in the vicinity of the control panel.



Version001_EN